

heydent

Experience the difference

2024

Dental-
Materials

Made in Germany





heydent
Experience the difference

Heydent GmbH
Viktor Frankl Str. 20, 86916 Kaufering
Germany
www.heydent.com



Composites

Page

Nano Ceramic Composite	HEYTEC-NC	3
Hybrid Composite	HEYTEC-BLEACH (X)	3
Nano Composite	HEYTEC-NR	4
Nano Hybrid Composite	HEYTEC-N	5
Micro Hybrid Composite	HEYTEC-MR	6
Universal Shade	HEYTEC-Universal	7
Bulk Fill Composite	HEYTEC-Bulk UD	8
Flow Composite Viscous	HEYTEC-Flow-V	9
Modeling Liquid	H-Moli	9

Adhesives

Bond, selfetching, 7th generation	H-ETCHBOND	10
Bond, 5th generation	H-BOND	10
Bond Activator	H-BOND Activator	10

Temporary Material

Temporary Crown and Bridge Material	HEYTEC Temp	11
Temporary Cement	HEYTEC Semi	13

Cement

Universal Resin Cement	HEYTEC CemX	12
Veneer Cement	HEYTEC Cem Veneer	12
Temporary Cement	HEYTEC Semi	13
Permanent Cement, selfetching	HEYTEC Cem SE	14

Core Build Up

Core Build Up Composite	HEYTEC Core-N	15
-------------------------	---------------	----

Bleaching

In Office Bleaching	JW NEXT Power Bleaching	16
Dam, liquid	JW Liquid Dam	16

Desensitizing

Silikon-Based Laquer	H-NANO Seal	17
----------------------	-------------	----

Impression Material


A-Silicon Soft	HEYTEC Soft,	18
A-Silicon SuHy (Super Hydrophilic)	HEYTEC Soft-H	18
A-Silicon Heavy	HEYTEC Heavy-X	19
A-Silicon Mono	HEYTEC Mono-X	19
A-Silicon Wash	HEYTEC Light-H	20
A-Silicon Bite Registration	HEYTEC Bite	20
A-Silicon Adhesive	HEYTEC Adhesive	20

HEYTEC-NC

HEYTEC Nano Ceramic is a ceramic-like, light-curing nano ceramic filling material and a multi-talent for long-lasting fillings, with significantly reduced sensitivity to ambient light and excellent physical values. **HEYTEC Nano Ceramic** is indicated for restorations of all cavity classes and for making inlays, onlays and veneers.

HEYTEC-NC

Nano Ceramic

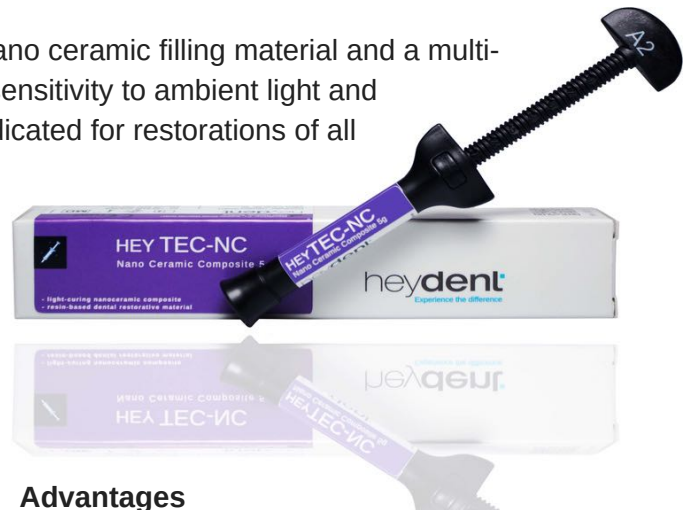

Order number	1x 5g syringe
HE011 - A1	HE017 - B2
HE012 - A2	HE021 - Inc
HE013 - A3	HE027 - BLE
HE014 - A3,5	HE028 - BLX
HE015 - OA2	HE073 - Set 8x5g
HE022 - OA3	HE078 - Set 8x5g+Bond
HE016 - B1	HE078 - Set 8x5g+Etchbond


Advantages

- abrasion resistance
- low shrinkage
- "Chameleon effect" for excellent aesthetics
- easy polishing
- long lasting shine
- comfortable modeling
- non sticky effect
- wide range of applications

TECHNICAL DATA

- Vickers Hardness 948 MPa
- Flexural Strength 157 MPa
- E-Module 16.575 MPa
- Compressive Strength 432 MPa
- Tensile Strength 56 Mpa
- Vol. shrinkage: 1,66%
- Radiopacity: 220%Al


HeyTEC-BLEACH-X BLEACH

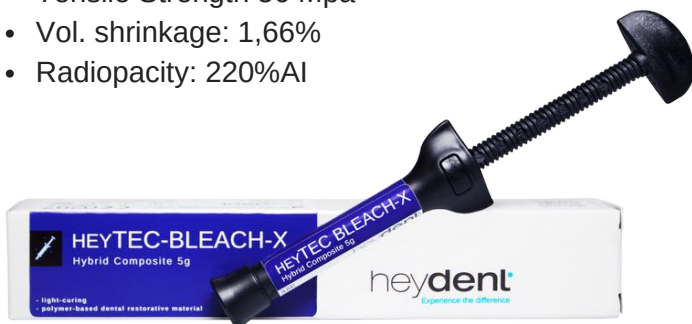

HEYTEC BLEACH & BLEACH-X is a ceramic hybrid composite material which is widely used for direct composite veneer cases. The excellent handling and material properties guarantee highly esthetic works. **HEYTEC BLEACH** is already white, but the extra bleach color is even whiter.



Order number	1x 5g syringe
HE025 - BLEACH	
HE026 - BLEACH-X	

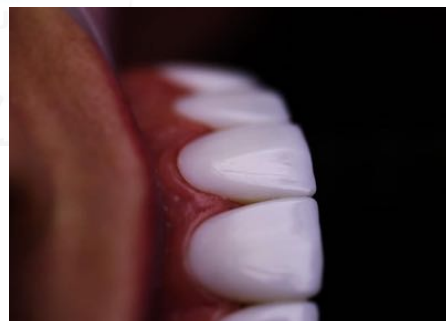

Direct Composite Veneers with HeyTEC BLEACH-X

HEYTEC-NR

HEYTEC NR is our new revolutionary Nano Composite in the Heydent restorative family.

HEYTEC-NR is a Multi-Talent for long-lasting fillings and with significantly reduced sensitivity to ambient light plus excellent physical values. Indicated both, for restorations of all cavity classes and for making inlays, onlays and veneers.



HEYTEC-NR

Nano Revolution

Compare to:

**Venus Diamond (Heraeus Kulzer),
IPS Empress Direct (Ivoclar Vivadent),
Filtek Supreme XTE (3M Espe),
Amaris (VOCO)**

Order number	1x 4g Syringe
HE511 - A1	HE521 - OA3
HE512 - A2	HE515 - B1
HE513 - A3	HE570 - Set 8x4g
HE514 - A3,5	HE578 - Set 8x4g+Bond
HE520 - OA2	HE076 - Set 8x4g+Etchbond

Technical Data

Vickers hardness	948 MPa
E-Module	16575 MPa
Compressive strength	432 MPa
Flexural strength	157 MPa
Tensile strength	56 MPa
Vol. shrinkage	1,66 %
Radiopacity	220%Al
Stability to ambient light	3 minutes



Features


Features	Benefits
Improved formula	Optimized working time
Nano technology	High physical values
Excellent mechanical properties	Longevity for long-lasting fillings
Wide range of shades	Natural looking restorations

HEYTEC-N

Nano Hybrid Composite **HEYTEC-N** is so convenient and easy to form, that its consistency is called „Wax-like“! **HEYTEC-N** enables the dentist to create highly aesthetic restorations with mechanical properties that guarantee the long-term reliable results.



HEYTEC-N Nano Revolution

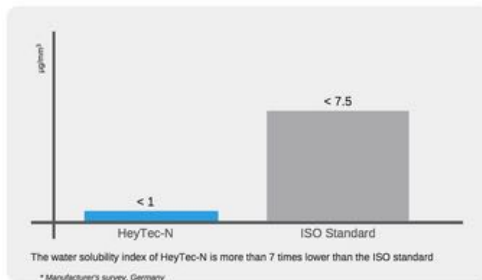


LONGTERM AESTHETIC FILLINGS

Thanks to its excellent resistance to compression, **HEYTEC-N** is able to withstand high loads on both occlusal and incisal surfaces! Furthermore, **HEYTEC-N** has an unusually low solubility in water (more than 7 times lower than the ISO standard)! **HEYTEC-N** is ideal for longterm aesthetic restorations.

ABRASIVE WEAR RESISTANCE

Studies on abrasive wear show that **HEYTEC-N** achieved one of the best results, among its competitors. After 200,000 cycles the surface of the filling shows only 38µm wear. **HEYTEC-N** restorations are extremely durable and safe for the opposite teeth.

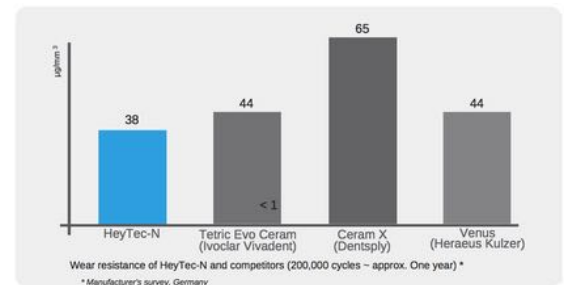

Order number	1x 4g Syringe
HE311 - A1	HE316 - B1
HE312 - A2	HE317 - B2
HE313 - A3	HE319 - D2
HE314 - A3,5	HE320 - D3
HE328 - A4	HE321 - Inc
HE315 - OA2	
HE322 - OA3	

Order number	20x 0,3g Compules
HE931 - A1	HE933 - A3
HE932 - A2	HE934 - A3,5



HEYTEC-MR

HEYTEC-MR is a universal light cure filling material with a high content of inorganic filler. **HEYTEC-MR** is perfect for direct restorations in both the anterior and posterior areas and for making inlays, onlays and veneers.

HEYTEC-MR is also suitable for core build-up and extended fissure sealing. The special combination of macro and micro fillers gives **HEYTEC-MR** its outstanding physical properties and makes it easy to polish. Thanks to these excellent qualities, this material is extremely plaque-resistant.

HEYTEC-MR is scientifically proven to be especially biocompatible and gentle to the gums.

TECHNICAL DATA

- Vickers Hardness 784 MPa
- Flexural Strength 157 MPa
- E-Module 11.173 MPa
- Compressive Strength 502 MPa
- Tensile Strength 50 Mpa

Order number	1x 4g syringe
HE411 - A1	HE427 - Incisial
HE412 - A2	HE425 - Bleach
HE413 - A3	HE426 - Bleach-X
HE414 - A3,5	HE470 - Set 8x4g
HE416 - B1	HE478 - Set 8x4g+Bond
HE417 - B2	HE476 - Set 8x4g+Etchbond

HEYTEC-MR

New microhybrid technology

ADVANTAGES

- Indicated for use in direct restorations of all cavity classes.
- Scientifically tested
- Highly advanced physical properties
- Excellent mechanical data
- Polishes well, homogenous surface, special particle structure
- Radiopaque
- Innovative colour scale, wide range of shades and highly- transparent incisal masses



HEYTEC-Universal HEYTEC-Universal BLOCK

HeyTEC Universal was developed for users who only want to use a shade-adaptive composite which can be used for all natural tooth shades for anterior and posterior restorations. **HeyTEC Universal** shows excellent physical properties and good handling features. The total filler content is 75 % by weight.

ADVANTAGES

- High transparency
- Radiopaque
- Excellent chameleon Effect
- One shade covers nearly all tooth colors





TECHNICAL DATA

- Vickers Hardness 790 MPa
- Flexural Strength 150 MPa
- E-Module 12.655 GPa
- Compressive Strength 479 MPa

Order number	1x 4g syringe
HE432 - Universal Color	
HE433 - Universal Block	

HEYTEC-Universal CHAMELEON effect

In case of heavily discolored cavities you can use **HEYTEC-Universal Block**, an opaque blocker, used for cavities in Class III or IV cases to cover the color and create a base for **HEYTEC-Universal**

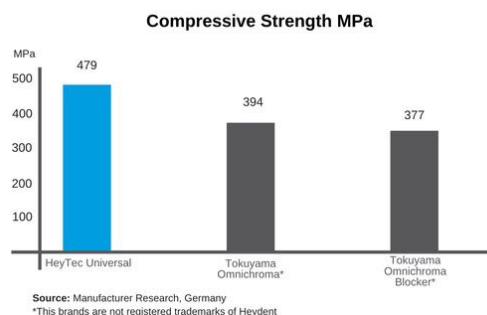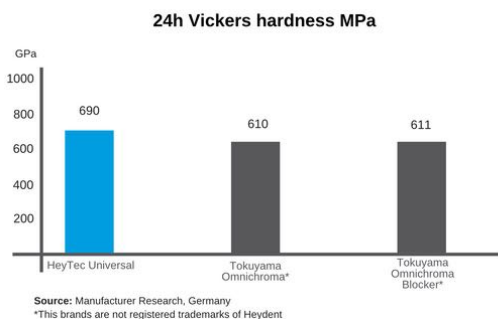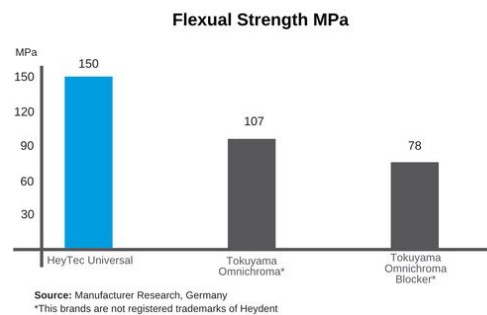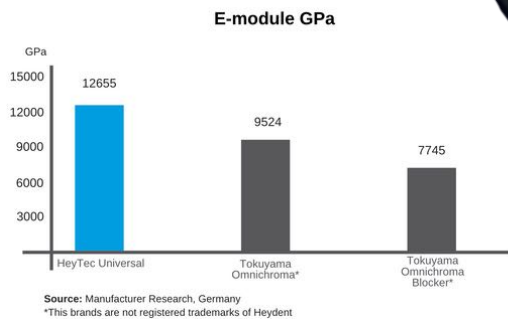

HEYTEC-Bulk

HEYTEC-Bulk is a flowable composite for direct fillings, using the Bulk Fill technique. Thanks to the optimal ratio of flowability to elasticity, it is quick and convenient to apply and shape.

HEYTEC-Bulk is intended for single-step reconstructions, but also layer technique.

The special monomer matrix ensures that the material adapts perfectly to cavity walls, even with the presence of undercuts.



CHAMELEON effect

Features


HEYTEC-Bulk

- up to 4 mm layers
- universal dentin shade with „chameleon” effect
- optimal ratio of flowability to elasticity
- excellent physical and mechanical properties
- high radiopacity (210% Al)
- Bis-Gma free


Benefits

HEYTEC-Bulk


- highly efficient handling
- saves time
- excellent adaption to cavity bottom and walls
- high quality, long-lasting posterior restorations
- perfect visual control
- optimised biocompatibility

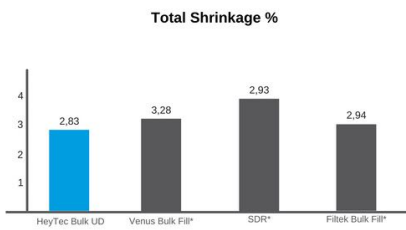

Source: Independent InnoMat GmbH, Tellow, Germany, 2012.
*This brands are not registered trademarks of Heydent.

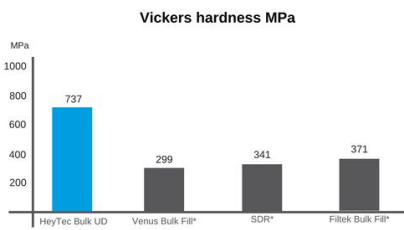


Source: Fraunhofer Institut, Dresden Germany 2012
*This brands are not registered trademarks of Heydent.

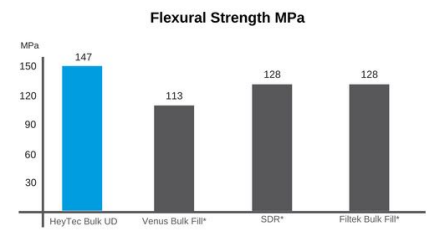


Source: Fraunhofer Institut, Dresden Germany 2012
*This brands are not registered trademarks of Heydent.

Order number	Universal Dentin
HE095 UD 2g syringe	
HE096 UD 10 x 0,2g compules	

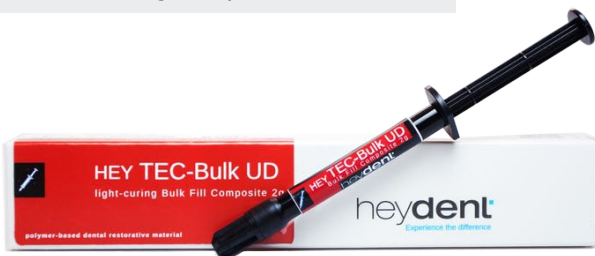

TECHNICAL DATA

- Bonding strength to dentine: 26 MPa
- Bonding strength to enamel: 24 MPa
- Vickers hardness 737 MPa
- Flexural strength 147 MPa
- E-Module 12.316 MPa
- Compressive strength 381 MPa

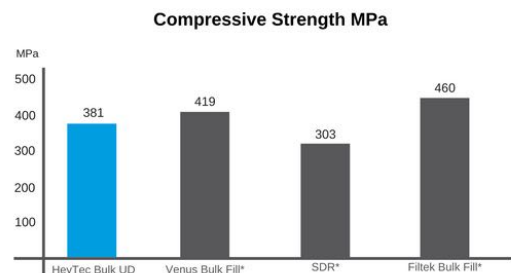

Source: Fraunhofer Institut, Dresden Germany 2012
*This brands are not registered trademarks of Heydent.

HEYTEC-Flow-V

HEYTEC-Flow-V is a new generation of a light-curing and radiopaque flow composite for direct restorations, relinings, as well as sealing fissures. **HEYTEC-Flow-V** is flowable as well as viscous when necessary and adapts nicely to cavity walls. The ideal flow characteristics simplify the restoration.

Features

HEYTEC - Flow-V

- Optimal flow behaviour
- excellent physical and mechanical properties
- high radiopacity (210% Al)
- Bis-Gma free

TECHNICAL DATA

- Vickers hardness 794 MPa
- Flexural strength 144 MPa
- E-Module 12.316 MPa
- Compressive strength 431 MPa
- Filler content approx. 77%
- Radiopacity 210% Al

Order number	1x 1,2ml syringe
HE611 A1	HE615 B1
HE612 A2	HE618 Clear
HE613 A3	HE619 Bleach
HE614 A3,5	

Benefits

HEYTEC - Flow-V

- Both flowable and viscous No running or dripping, easy and highly precise application. Even penetrates into areas that are difficult to reach
- Stable and long-lasting restorations
- perfect visual control
- optimized biocompatibility



HEYTEC-MoLi



Modeling Liquid

HEYTEC-MoLi is the solution to build up resin composite (RC) restorations, easy, fast & efficient.

Order number	1x bottle a 5ml
HE004	

H-ETCHBOND


H-ETCHBOND is a selfetching, light-curing one-component, single-layer bonding system of the seventh generation! Developed to simplify the bonding procedure, while ensuring excellent quality and strong bonding.

7 Generation


7.Generation

One of the strongest bonds - test and compare!

Adhesion strength on enamel

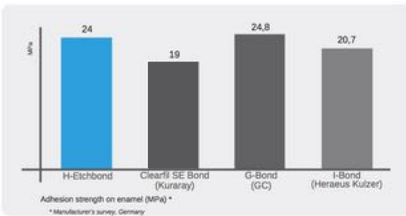

Adhesion strength on dentin

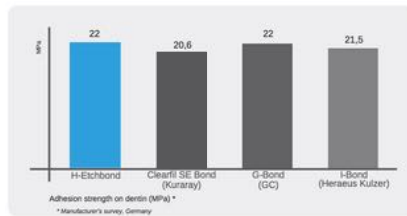

Amongst the most popular selfetching systems, **H-ETCHBOND** is characterized by its outstanding adhesion properties. Only one layer of **H-ETCHBOND** is required to obtain a strong and long-lasting bond of light-curing composites to enamel and dentin.

Order number 1x bottle a 5ml
HE003

H-BOND

H-BOND is a light-curing system of the Prime / Bond type of the 5th generation and offers a strong bond of all composite materials to dentin and enamel! **H-BOND** has been developed for light-curing materials, but when combined with **H-BOND Activator**, **H-BOND** becomes a chemically or dual-curing system that is perfect for root canals! It can also be used with chemically and doubly hardening materials.

H-BOND Activator



Order number 1x bottle a 5ml
HE001

Order number 1x bottle a 3ml
HE002

HEYTEC Temp

Highly esthetic 4:1 composite material for provisional crowns and bridges inlays, onlays and veneers

Properties

- High fracture resistance also in very thin layer such as edges of crowns or thin veneers
- Very aesthetic appearance with high colour stability and natural fluorescence
- Material of the latest generation with mixing ratio 4:1
- Very high mechanical durability and excellent elasticity
- Optimal for long bridges with high durability
- Excellent fit due to minimal polymerisation shrinkage
- High abrasion resistance
- Excellent polishability
- Safe and quick handling with self mixing cartridge system
- Direct chairside application due to low curing temperature
- Final preparation after only 5 min

TECHNICAL DATA

- Working time at 23 ° C / 73 ° F 0:50 min
- Shore-D-hardness 75
- Barcol-hardness 36
- Flexural strength 88 MPa
- Compressive strength 304 MPa
- Diametral tensile strength 44 MPa
- Elasticity Modulus 2,7 GPa

Temporary Crown and Bridge Material

Easy application and preparation

HEYTEC Temp can be mixed and applied automatically, which ensures the homogenous consistency of the material. Formation of air bubbles is avoided and the risk of cracks is minimized. Controlled chemical curing of the material allows the convenient application

inside and outside of the patient's mouth.

The easy excess removal with a scalpel or scissors is done, before the final preparation with the rotary instrument, during the elastic phase.

Order number	1x 50ml, 10x MixTips
HE089 - A1	
HE056 - A2	
HE090 - A3	
HE093 - Bleach	



HeyTEC CemX

HeyTEC CemX is a dual-cure, highly radio-opaque, fluorescent and highly abrasion-resistant flowable self-adhesive resin cement for various types of indications

Universal Resin Cement

HeyTEC CemX is for various indications such as all-ceramic, composite or metal inlays, onlays, crowns and bridges, posts, 3-unit inlay/onlay bridges, all-ceramic, 2- or 3-unit maryland bridges and veneers plus **HeyTEC CemX can also used for core build up.**

HeyTEC CemX is the latest generation of self-adhesive resin cements.

TECHNICAL DATA

- Vickers hardness - >549 MPa
- Compressive strength - >290 MPa
- Flexural strength - >118 MPa
- E-Module - >12039 MPa
- Tensile strength - 45 MPa
- Radiopacity - 210%
- Polymerisation - 7-9 Minutes



Order number	1x 8g syringe
HE033	semi transparent
HE034	A1
HE035	A2
HE036	A3
HE037	B1

Features

HeyTEC - CemX

- One material for three indications
- High viscosity
- Automix cartridge
- Highly radiopaque
- High filler content
- Excellent physical data
- Dual-curing
- Bis-Gma free

Benefits

HeyTEC - CemX

- Optimal "interface" between pin-cement-dentin-crown
- Highly abrasion resistant
- Very comfortable handling
- Long-lasting restorations and extremely low shrinkage
- Direct and ergonomic application/no mixing necessary
- Curing can be accelerated with polymerisation light
- X-ray check possible
- Optimised biocompatibility

HEYTEC Cem Veneer

HEYTEC Cem Veneer is a light-cure, single component, highly radiopaque luting material for cementing porcelain, composite and ceramic veneers and for cavity lining.

HEYTEC Cem Veneer is BisGMA, Bis EMA and TEGDMA free and highly viscous. HEYTEC Cem Veneer offers easy handling and excellent color stability for highly esthetic veneer cases.

- Filler Content 80%
- Compressive strength - >362 MPa
- Flexural strength - >140 MPa
- Film thickness - 28um
- Viscosity - 590 Pas
- Radiopacity - 210%
- E-Module - >14900 MPa



Order number	1x 2g syringe
HE044	A2
HE041	BLX

HEYTEC Semi

HEYTEC Semi is a dual-curing, semi-permanent cement for the cementation of crowns and bridges on implants and natural teeth. No chemical adhesion [0 MPa] is compensated by controlled expansion of the cement. This ensures optimal fastening, and removal of the restoration without any damage.

Temporary Cement

TECHNICAL DATA

- Compressive strength - 40 MPa
- Film thickness - 20 µm
- Processing time (23° C / 73° F) - 2:00 min.
- Setting time (self-curing) - approx. 4:00 min.
- Light curing - 40 sec per area



Order number 1x 5ml, 5x mix tips
HE032


Excellent "snap-off" for revisions

Revision or prosthetic reconstruction removal is safe and easy. **HEYTEC Semi** application does not require etching, conditioning or bonding, which eliminates the chemical connection between the abutment and crown. Cement residue can be easily removed and the construction recemented if necessary.

Biocompatible and antibacterial

HEYTEC Semi is eugenol free, has a neutral taste and smell for patient comfort and added zinc oxide causes an antibacterial effect.

Compressive Strength MPa

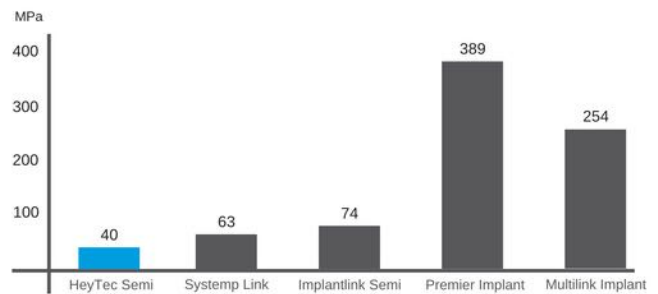


Source: Fraunhofer Institut, Dresden Germany
*This brands are not registered trademarks of Heydent

Thin layer and controlled expansion

precise fit is guaranteed by a thin layer of **HEYTEC Semi** (<20). Controlled volume expansion of the cement (eliminates the shrinkage). All this gives an ideal marginal sealing.

Expansion of HeyTec Semi vs. shrinkage of other cements

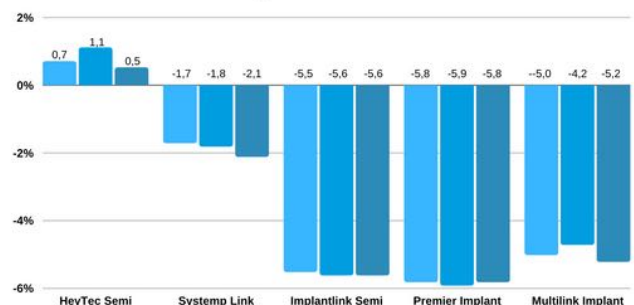

Source: Manufacturer Research, Germany

Flexible structure for the bone and implant protection

Stress-absorbing properties are much higher due to the retained flexibility of **HEYTEC Semi** (40MPa). This creates a much better environment for tissue regeneration.

HEYTEC Cem SE

HEYTEC Cem SE is a selfadhesive, dualcure composite material for final cementation of crowns, bridges, inlays/onlays restorations and posts. **HEYTEC Cem SE** provides high strength adhesion to dentin, enamel, zirconium (16MPa!), ceramics, composites, titanium and metal.

Selfetching Cement

Strong adhesion without etching, conditioning or bonding

The entire procedure of cementation requires just a few steps, with no risk of overetching and sensitivity. **HEYTEC Cem SE** is self neutralizing with an initial low ph, which increases after application to ph 5.5 after some hours. The rise of ph prevents post cementation sensitivity.



Order number 1x 5ml, 10x mix tips
HE031


Perfect fit reconstruction and predictable work

A thin film - 12µm - provides precise cementation while delicate expansion of the cement eliminates its shrinkage from the margins. This ensures excellent marginal integrity. The dual cure system allows pre-curing of the positioned restoration, as well as easy removal of excess in its initial flexible phase.

Excellent color stability and natural fluorescence

Proven color stability is satisfying even for transparent ceramics, this guarantees a long-lasting aesthetic effect. An additional advantage is the significant release of fluoride, which results in strengthening the surrounding tissues.

HeyTec Cem SE
Self neutralization after 12h

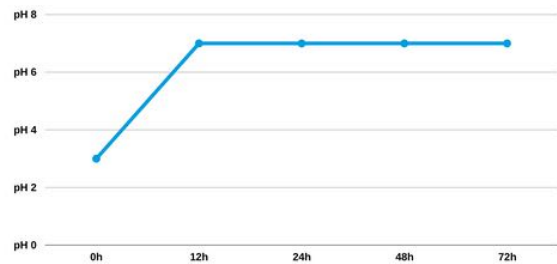


Source: Manufacturer Research, Germany

TECHNICAL DATA

- Shear bond strength on enamel 9 MPa / dentin 12 MPa
- Compressive strength 180 MPa
- Setting time (self cure) 4:00 min.
- Light-curing (optional) 20 sec.

Film thickness comparison µm

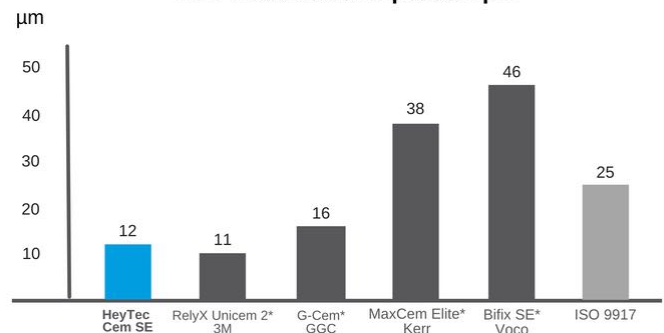

Source: Manufacturer Research, Germany
*This brands are not registered trademarks of Heydent

HEYTEC Core-N

HEYTEC Core-N is a dual-cure flowable composite material with nano particles (nano-zirconium dioxide, nano-calciumfluoride) for core-build-up and root post cementation.

Reliable and durable **Core Build Up**

NaNo-technology

The exceptional mechanical properties of HEYTEC Core-N is caused by the innovative molecule structure. The material is highly resistant against compression and bending. These features make HEYTEC Core-N the perfect replacement of the natural dentin.

Secure and predictability treatment



HEYTEC Core-N can be shaped exactly like dentin, which makes the work easy and predictable. Extra safety, is guaranteed by no temperature increase during material polymerization, especially when reconstruction vital teeth

Dentin remineralization and x-ray visibility

The content of nano-particles of calcium fluoride allows dentin remineralization, and nano-zirconia particles generate excellent x-ray visibility.



HEYTEC Core-N is a multi purpose material for core build-up and for the cementation of endo posts. Together with a root post, it creates a homogeneous and durable structure for the permanent tooth reconstruction.

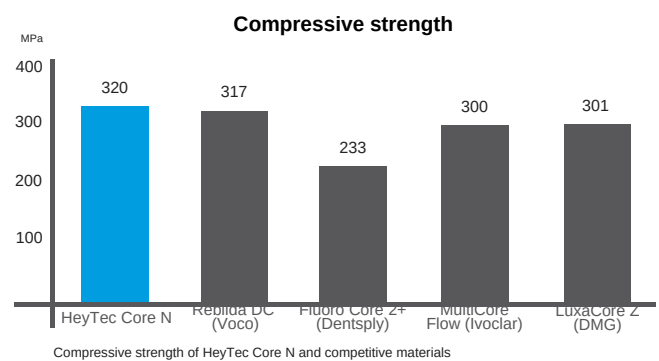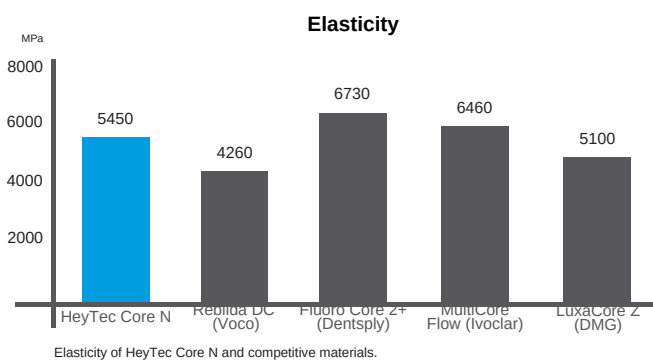

TECHNICAL DATA

- Bending strength (dual-hardening mode) - 145 MPa
- Flexural strength (self-hardening mode) - 129 MPa
- Compressive strength - 320 MPa
- Elasticity (dual-hardening mode) - 7.4 GPa
- Elasticity (self-hardening mode) - 5.4 GPa

Order number 1x 5ml, 5x mix tips
HE049

JW NEXT Power Bleaching

JW NEXT Power Bleaching is a two-component chair-side bleaching gel in a double-bar cartridge with 35 % content of active ingredient (hydrogen peroxide) for quick bleaching and brightening of discolored and devitalized teeth.

JW NEXT Power Bleaching can be used without any light source, but gets better results with lamps and best results with a diode-laser!



Order number	1x SET*
1020JW	

* 1x 5g Double cartridge, 1x 1,5g JW Liquid Dam, 3x Mixing Tips, 1x Liquid Dam tip



Clinical case from Dr Eugenia Anagnostaki

Features

JW NEXT Power Bleaching

- High content of active ingredient
- Ideal consistency
- Automix cartridge
- Neutral pH-value
- Red color

Benefits

JW NEXT Power Bleaching

- Quick bleaching and brightening
- Long-lasting result
- Easy application
- No running, good contrast on the tooth
- Consistent quality
- No mixing necessary

JW Liquid Dam



JW Liquid Dam is a one-component material to cover the gingiva during treatment with highly active in-office bleaching gels. **JW Liquid Dam** provides optimal viscosity, good flow behavior, very low heat development and only 20 sec polymerization time. It's color is blue for a very good contrast to the gingiva.

Order number	2x 1,5g syringes
1003	

H-NANO Seal

H-NANO Seal is an innovative protective lacquer, based on a patent-protected silicone polymer!
H-NANO Seal forms a very strong adhesion on the enamel or dentin. It is not necessary to etch before.
Thanks to the innovative formula, the material stays on the tooth surface for about 12 months.

Desensitizer



12 months

protection up to 12 month

Features

H-NANO Seal

no need to polimerize

no etching

Silikonpolymer

- Ensures excellent adhesion to enamel and dentin
- Transparent - does not cause any stains
- Can be used on bleached teeth

NANO Fluorapatit

- Fills enamel microfractures
- Closes open tubules

Aminfluorid Olaflur

- Immediately releases large amounts of fluoride
- Accelerates the process of remineralization

Benefits

H-NANO Seal

- caries prophylaxis for children and adults, as well as patients with braces
- treatment of teeth after a professional cleaning
- tooth protection after teeth whitening
- reduction of hypersensitivity caused by open tubules
- remineralization of hard tooth tissues and closing of microfractures in the enamel
- preliminary sealing of fillings and restorations

Order number

1x set*

HE008

* 1x 5ml bottle, 20x disposable brushed, applicaton trays

HEYTEC Soft

Kneadable A-Silicone for the double mixing technique. **HEYTEC Soft** has a particularly soft mixing consistency. The impression is taken without excessive compression pressure and can be made in both prefabricated and individual impression trays.

- Consistency: putty
- Color yellow
- Linear shrinkage: <0.1%
- Mixing time: 45 sec.
- Processing time: 2:20 minutes
- Curing time (in the mouth): <3:00 minutes
- Final hardness: approx. 65 shore A
- Standard: ISO 4823, type 0



Order number

HD060 2x 300 ml Jar A+B, 2 x Spoon

HEYTEC Soft-H

Kneadable A-silicone for the double mixing or double impression technique. Due to the smooth initial consistency and final hardness of 67 shore A, **HEYTEC Soft-H** (superhydrophil) is universally applicable. The extremely strong hydrophilicity develops a pronounced affinity to all wetted areas and guarantees the highest possible surface quality of the impression.

- Consistency: soft putty
- Color: gray - green
- Linear shrinkage: <0.1%
- Mixing time: 45 sec.
- Processing time: 2:20 minutes
- Curing time (in the mouth): <3:00 minutes
- Final hardness: approx. 67 shore A
- Standard: ISO 4823, type 0



Order number

HD068 2x 300 ml Jar A+B, 2 x Spoon

HEYTEC Mono-X

Monophase Precision A-Silicone with excellent stability and exceptional image sharpness.

Areas of application: Inlay-Onlay crowns and bridges, fixation impressions, telescopic and attachment work up to complex implant impressions.

- Consistency: heavy
- Color: aubergine
- Linear shrinkage: <0.1%
- Mixing time: not applicable
- Processing time: 2:00min
- Curing time (in mouth): <3:00 min
- Finalhardness: approx. 70 shore A
- Standard: ISO 4823, Type 1

Order number

HD064 1x 380 ml 5:1 Cartridge, 10 Mix-Tips



HEYTEC Heavy-X

Precision impression material for machine, dynamic mixing. **HEYTEC Heavy-X** is suitable for all precision impression techniques in which outstanding hydrophilic properties, excellent dimensional stability, high tensile strength and high resilience are important.

Areas of application: Pre-impression material for all two-phase impression techniques.

- Consistency: heavy
- Color yellow
- Linear shrinkage: <0.1%
- Mixing time: not applicable
- Processing time: 1:30 min
- Curing time (in mouth): <2:30 min
- Final hardness: approx. 65 shore A
- Standard: ISO 4823, Type 1

Order number

HD065 1x 380 ml 5:1 Cartridge, 10 Mix-Tips



HEYTEC Light-H

Addition-crosslinking, extremely fine-flowing correction impression material with superhydrophilic properties for all impression techniques. With the avoidance of air inclusions, it flows to all critical points and also reproduces the finest contours. Residual liquids are washed out by the compression pressure and the superhydrophilic properties. The thixotropic property prevents inadvertent flow of the already applied material from the preparation.

- Consistency: putty
- Color yellow
- Linear shrinkage: <0.1%
- Mixing time: 45 sec.
- Processing time: 2:20 min
- Curing time (in mouth): <3:00 min
- Final hardness: approx. 65 shore A
- Standard: ISO 4823, type 0



Order number

HD061 2x 50 ml Cartridge, 12 Mix-Tips

HEYTEC Bite

Addition-crosslinking, superhard, fast-setting silicone material for bite registration. With a final hardness of approx. 90 shore A, the material is not only suitable for the registration of full or partial toothed jaws, but also for bite registration in the total prosthesis. The post-processing with knife, stone or milling machine is possible without any problem.

- Consistency: heavy
- Color: magenta
- Linear shrinkage: <0.1%
- Mixing time: not applicable
- Processing time: 15 seconds
- Curing time (in mouth): 45 seconds
- Final hardness: approx. 90 shore A
- Standard: ISO4823, type 0



Order number

HD080 2x 50 ml Cartridge, 12 Mix-Tips



HEYTEC Adhesive

Liquid spoon adhesive for all commercially available A-silicones.



Order number 1x 10ml

HD062

 **WhatsApp**
Click or Scan
to chat with us! 



heydent

Experience the difference



Heydent GmbH, Viktor Frankl Str. 20,
86916 Kaufering, Germany

www.heydent.com

info@heydent.com

Tel: +49 8191 9737804 Fax: +49 8191 9855137